 **BD Veritor™ Plus Analyzer**
**Instructions for Use/Instructions d'utilisation/
Gebrauchsanweisung/Istruzioni per l'uso/Instrucciones de uso/Bruksanvisning**

  **R_x Only**

 Becton, Dickinson and Company
7 Loveton Circle
Sparks, Maryland 21152 USA

 Benex Limited
Pottery Road, Dun Laoghaire
Co. Dublin, Ireland


Australian and New Zealand Sponsors

Becton Dickinson Pty Ltd.
66 Waterloo Road
Macquarie Park NSW 2113, Australia

Becton Dickinson Limited
14B George Bourke Drive
Mt. Wellington Auckland 1060, New Zealand

 bd.com/e-labeling  256066
8091498(08) 2021-02

Affix the USB Port Unlock label here.

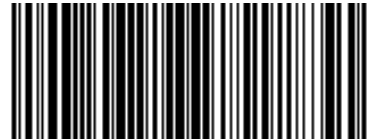8091498(08)

ATTENTION: To use this BD Veritor™ Plus Analyzer with BD Veritor™ System assays, you must download the most current package inserts and Quick Reference Instructions from bd.com/e-labeling.

Change History/Historique des modifications/ Bisherige Änderungen/Cronologia delle modifiche/ Historial de cambios/Revideringshistorik

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LIMITED WARRANTY

This warranty gives you specific legal rights. Additionally, you may have other rights that vary from region to region.

The BD Veritor™ Plus Analyzer, BD Veritor™ InfoScan Module, BD Veritor™ Plus Connect, and BD Veritor™ InfoWiFi module are warranted to be free from defects in materials and workmanship for a period of two years following delivery. BD's sole responsibility under this warranty shall be to repair or replace any instrument or its components (except for expendable supplies) which under normal operating conditions, prove to be defective within two years of delivery.

BD will furnish new or re-manufactured components upon its option. All replacements shall meet new part specifications. Replaced components become the property of BD.

It is understood that the equipment covered by this warranty has been installed in accordance with the recommendations and instructions in the *BD Veritor™ Plus Analyzer Instructions for Use Manual*.

Any damage to a BD Veritor™ Plus Analyzer resulting from the failure of the owner to maintain reasonable care and precautions in the operation and maintenance of the system will void this warranty and terminate the obligations of the manufacturer as stated herein.

This warranty is in lieu of all other warranties, whether express or implied, including but not limited to, warranties of merchantability, or fitness for a particular use. In no event will BD Diagnostic Systems be liable for indirect, incidental or consequential damages.

No part of this publication may be reproduced, transmitted, transcribed, stored in retrieval systems, or translated into any language or computer language, in any form or by any means, electronic, mechanical, magnetic, optical, chemical, manual, or otherwise, without the prior written permission of BD, 7 Loveton Circle, Sparks, Maryland, 21152, USA.

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BD, the BD Logo, Synapsys, and Veritor are trademarks of Becton, Dickinson and Company or its affiliates. All other trademarks are the property of their respective owners. © 2021 BD. All rights reserved.

Table of Contents

1. Introduction	
1.1 Intended Use	9
1.2 Summary and Explanation	9
1.3 Principles of the Procedure	10
1.4 Use of this Manual	11
1.5 Conventions	11
1.5.1 User Interface	11
1.5.2 Notes, Cautions, and Warnings	11
1.5.3 Summary of Cautions and Warnings	11
2. Installation	
2.1 General	12
2.2 Analyzer Specifications	12
2.3 Analyzer Setup	13
2.4 Printer Setup – Optional	15
2.5 Module and Accessory Specifications – Optional	15
2.6 BD Veritor InfoScan Module Setup – Optional	16
2.7 BD Veritor Plus Connect Set Up – Optional	17
2.8 BD Veritor InfoWiFi Module Setup – Optional	17
3. Controls and Indicators	
3.1 General	18
3.1.1 Coaxial Power Port	18
3.1.2 Power Button	18
3.1.3 System Indicators	19
3.1.4 Test Device Slot	20
3.1.5 USB Port	20
4. Operation	
4.1 General	20
4.2 Using the BD Veritor Plus Analyzer – Analyze Now Operation	20
4.3 Using the BD Veritor Plus Analyzer – Walk Away Operation	22
4.4 Kit Quality Control (QC) Test Procedure – Analyze Now Operation	23
4.5 Kit Quality Control (QC) Test Procedure – Walk Away Operation	25
4.6 Workflow Configuration Procedure – Optional	27
4.7 Workflow Configuration Barcodes	27
4.8 Data Download to a PC	28
5. Maintenance	
5.1 General Procedures for Maintenance and Servicing	29
5.1.1 BD Veritor System Verification Cartridge	29
5.1.2 Cleaning	30
5.2 Servicing	30
5.3 Transportation and Storage	30
6. Troubleshooting	
6.1 Instrument Service	30
6.2 Forced Power Down	30
6.3 Troubleshooting Aids	30
6.4 Result Messages	31

7. Configuration Barcodes 38

Translations/Traductions/Übersetzungen/Traduzioni/Traducciones/Översättningar

Français 41
Deutsch 75
Italiano 109
Español 145
Svenska 179

International Contacts/Contacts internationaux/Internationale Ansprechpartner/Centri di assistenza internazionali/Contactos internacionales/Internationella adresser

International Contacts/Contacts internationaux/Internationale Ansprechpartner/Centri di assistenza internazionali/Contactos internacionales/Internationella adresser 213

Symbol Glossary/Glossaire des symboles/Symbolverzeichnis/Legenda dei simboli/Glosario de símbolos/Symbolförklaring

Symbol Glossary/Glossaire des symboles/Symbolverzeichnis/Legenda dei simboli/Glosario de símbolos/Symbolförklaring 215

1 Introduction

1.1 Intended Use

The BD Veritor™ Plus Analyzer is for use with BD Veritor™ System Test Devices.

1.2 Summary and Explanation

The BD Veritor™ Plus Analyzer is intended to provide rapid test results in near-patient settings.

The BD Veritor™ Plus Analyzer is a digital immunoassay instrument that is a component of the BD Veritor™ Plus System. The Analyzer supports the use of different assays by reading an assay-specific barcode on the test device. Depending on the configuration chosen by the operator, the instrument communicates status and results to the operator via a liquid crystal display (LCD) on the instrument, a connected printer, or through a secure connection to the facility's information system.

A BD Veritor™ Plus Analyzer can be identified by the image below:

The BD Veritor™ Plus System features a modular architecture that allows users to customize the system to meet their needs. Those who wish to include data capture and/or secure transmission features may expand the functionality of the system using optional accessories as described below. Please speak to your BD representative about the availability of each option, as all products may not be available in all markets.


BD Veritor™ InfoScan Module – Optional

The BD Veritor™ InfoScan module is inserted into the BD Veritor™ Plus Analyzer to add the capability of reading specimen identification, operator identification, reagent lot information, reagent expiration date, and modifying the on-screen display language of the BD Veritor™ Plus Analyzer. With the BD Veritor™ InfoScan module, users can unlock the unit's data storage to download test information to a connected computer over a USB connection. A BD Veritor™ Plus Analyzer configured with a BD Veritor™ InfoScan module can be identified by the label on its underside, shown in the image below.

BD Veritor™ Plus Connect – Optional

When implemented with BD Synapsys™ Informatics Solution, BD Veritor™ Plus Connect enables BD Veritor™ Plus Analyzers equipped with a BD Veritor™ barcode scanning module to securely transmit test results to a facility's information system via an Ethernet connection. The Connect software runs on a separate dedicated computer networked to the BD Veritor™ Plus Analyzer with the USB cable provided with the BD Veritor™ barcode scanning module. This is shown schematically below.

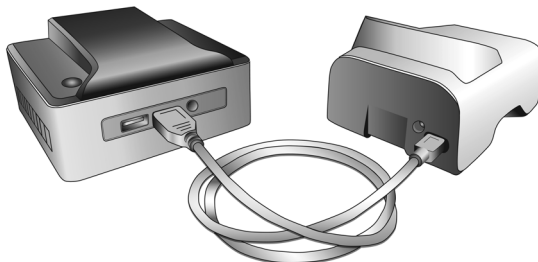BD Veritor™ InfoWiFi Module - Optional

The BD Veritor™ InfoWiFi module provides the same functional features as the BD Veritor™ InfoScan module and adds wireless communication capability. Once configured, the BD Veritor™ InfoWiFi module can connect to a local wireless network and transmit test results to the BD Synapsys™ Informatics Solution. An additional BD Synapsys™ Informatics subscription is required.

A BD Veritor™ Plus Analyzer configured with a BD Veritor™ InfoWiFi module can be identified by the label on its underside, shown in the image below.

1.3 Principles of the Procedure

The BD Veritor™ Plus Analyzer is an optoelectronic instrument that uses image processing and a specific algorithm to determine the presence or absence of any target analyte(s).

The BD Veritor™ Plus Analyzer detects the insertion of a BD Veritor™ Test Device and reads the device barcode to detect the assay type. Once it registers the test type, the Analyzer uses a reflectance-based measurement method to evaluate the line signal intensities on the assay test strip. The Analyzer determines the line intensity at each of the spatially-defined test and control line positions, interprets the results using a scoring algorithm, and reports a positive, negative, or invalid result based on pre-set thresholds. A liquid crystal display (LCD) on the instrument communicates status and results to the operator.

Because the BD Veritor™ Plus Analyzer and Test Devices are designed as an integrated™ system, it is not possible to replicate the result determination visually. Therefore assay result determination by operator visual readings of the BD Veritor™ System Test Devices is not recommended.

The instrument is maintenance-free and is configured with a single button design (power switch) for ease of operation.

1.4 Use of this Manual

This manual is designed as an integral part of instrument operation for trained laboratory personnel, technologists, supervisors, and other personnel who operate the BD Veritor™ Plus Analyzer instrument on a regular basis. Every attempt has been made to include all information that would be required during normal use and maintenance of the system. Should a question arise which is not answered in this manual, please contact the following parties:

In North America:

Technical Services 1.800.638.8663

International contacts are listed on page 213, or contact your local BD representative.

Comments or recommendations on this manual may be expressed via e-mail to LabelingDD@bd.com.

1.5 Conventions

1.5.1 User Interface

System prompts and messages are shown in CAPITAL LETTERS (e.g., INSERT TEST DEVICE).

1.5.2 Notes, Cautions, and Warnings

Throughout this manual, important information is presented in boxes offset from the regular text, and is labeled as either a NOTE, CAUTION, or WARNING. These messages are formatted as shown below and bear the following significance:

NOTE

Important information about instrument use worthy of special attention is presented as a NOTE.

CAUTION

Information on an activity which could potentially cause damage to the instrument is presented as a CAUTION.

WARNING

INFORMATION ON AN ACTIVITY WHICH COULD POTENTIALLY CAUSE INJURY TO THE USER IS PRESENTED AS A WARNING.

1.5.3 Summary of Cautions and Warnings

- For in vitro Diagnostic Use.
- Do not use after expiration date.
- Dispose of all BD Veritor™ kit components (test devices, reagents and control swabs) as biohazardous waste in accordance with local requirements after use.
- Dispose of expired BD Veritor™ Plus Analyzers as biohazardous waste and according to local environmental electronic waste disposal regulation.
- In the EU, after cleaning expired BD Veritor™ Plus Analyzers (see information in section 5) and placing the instrument in quarantine for seven days, dispose of them according to local environmental electronic waste disposal regulations.
- Do not eat, drink, or smoke in the area where the specimens and kits are handled.
- Keep the BD Veritor™ Plus Analyzer, BD Veritor™ System Test Devices and BD Veritor™ Verification Cartridge out of reach of children.
- Wear protective clothing such as laboratory coats, disposable gloves, and eye protection when specimens are processed.
- Humidity and temperature can adversely affect results; therefore, operate the BD Veritor™ Plus Analyzer indoors at room temperature (15–30 °C).
- Ensure that the BD Veritor™ Plus Analyzer is not in direct sunlight or exposed to a bright light.
- The BD Veritor™ Plus Analyzer is designed for use with BD Veritor™ System Test Devices. Do not insert any foreign materials into the BD Veritor™ Plus Analyzer.
- The BD Veritor™ Plus Analyzer may contain sensitive patient information and should not be stored in a publicly accessible location.
- BD recommends against using patient identifiers that can be readily traced back to a specific patient, such as a patient name or social security number.

- BD recommends affixing the unlock barcode provided with the BD Veritor™ Plus Analyzer to this manual and storing the manual in a secure location.
- Disposable materials and reagents should only be used within their expiration period. BD recommends that expired materials never be used.
- The BD Veritor™ Plus Analyzer contains lithium coin cell and lithium-ion cell batteries. Do not incinerate. Follow local disposal regulations. Separate collection for waste of electrical and electronic equipment. Do not dispose of battery in municipal waste. The symbol indicates separate collection for battery is required. Dispose of or recycle in accordance with national and local requirements.
- AC ADAPTER WARNING: RISK OF FIRE! SHOCK HAZARD!: To prevent electric shock, do not remove the cover. No user-serviceable parts inside. This unit contains hazardous voltages and should only be opened by a trained and qualified technician. Do not use if bent, cracked, modified or compromised. Use only the AC power adapter provided with the Analyzer.

Pathogenic microorganisms, including hepatitis viruses and Human Immunodeficiency Virus, may be present in clinical specimens. Standard Precautions and institutional guidelines for the handling of biohazardous materials should be followed in handling all items contaminated with blood and other body fluids.

Protection provided by this equipment may be impaired if the equipment is used in a manner not consistent with the instructions in this manual.

- EU Only: users should report any serious incident related to the device to the Manufacturer and National Competent Authority. Outside EU: Contact your local BD representative for any incident or inquiry related to this device.

2 Installation

2.1 General

This section provides specifications for installation and setup of the BD Veritor™ Plus Analyzer, plus the optional BD Veritor™ barcode scanning modules as well as the BD Veritor™ Plus Connect accessory. The following major topics are included:

- Analyzer Specifications
- Analyzer Setup
- Printer Setup
- Module Specifications
- Module Setup

WARNING

PROTECTION PROVIDED BY THIS EQUIPMENT MAY BE IMPAIRED IF THE EQUIPMENT IS USED IN A MANNER NOT CONSISTENT WITH THE INSTRUCTIONS IN THIS MANUAL.

2.2 Analyzer Specifications

Physical Specifications

Length	14.3 cm
Width	9.0 cm
Height (including feet)	7.6 cm
Weight (Analyzer)	300 g
Weight (Analyzer and installed module)	350 g

Electrical Specifications	
Supply Voltage	<p>The BD Veritor™ Plus Analyzer is provided with an auto-switching power supply with recommended input of 100–240 VAC at 50–60 Hz, 0.6 A.</p> <p> The Analyzer power requirement is 5.0 VDC ± 5%, 1,500 mA.</p>

Environmental Specifications	
Operating Conditions	
Ambient Temperature	15.0–30.0 °C (59–86 °F)
Ambient Humidity	Up to 85% RH, non-condensing
Ambient Light	0–3,000 lux uniform illumination
Ambient Environment	Indoor Use Only
Other	
Overvoltage Category	II (per IEC 60664)

Other Specifications	
Display Languages (BD Veritor™ barcode reader equipped module required for configuration)	English, French, German, Italian, Spanish, Swedish
Lifetime Number of Tests Performed	10,000 maximum*
Lifetime Number of Months from Date of First Use	24 months
Lifetime Maximum Shelf-Life from Date of Manufacture (regardless of use)	34 months
Operation Modes	Analyze Now mode Walk Away mode
* Effective in Version 5.50 firmware and later	

2.3 Analyzer Setup

Box contents


- 1 BD Veritor™ Plus Analyzer
- 1 Instructions for Use Manual
- 1 USB Port Unlock Barcode adhesive label (place label in a safe place)
- 1 compact AC power adapter with region specific blade adapters
- 1 BD Veritor™ Verification Cartridge

NOTE

If any items are missing, contact your BD representative. (See page 213 for contact information.)

Analyzer Unpacking and Setup

- 1 Before using the BD Veritor™ Plus Analyzer, visually inspect the contents for damage. If damage is apparent, contact your BD representative or BD Technical Services.
- 2 The BD Veritor™ Plus Analyzer ships with an adhesive-backed USB Port Unlock Barcode label, which is used on Analyzers that are configured with a BD Veritor™ barcode scanning module. BD recommends affixing this label to the space indicated on the cover page of this manual or storage in a secure location of your choice for future use. Instructions for using the USB Unlock Barcode can be found in section 4.6. The optional modules may not be available in all regions; please contact your local BD representative for availability.
- 3 Remove the Analyzer and power adapter from the protective wrapping. BD recommends retaining the packing materials in a safe location. The provided power adapter rated output is 5V 3A, 15W (Refer to the figure below).
- 4 Refer to the figure below. The North American plug is used to illustrate the installation of region specific blade adapters onto the power supply. Select the blade adapter that is appropriate for your region.
- 5 Slide the retaining clip (1) down on the power adapter.
- 6 Insert the notched end (2) of the blade adapter into the top of the power adapter at a 30–60 degree angle.
- 7 Push the bottom (3) of the blade adapter into the power adapter.
- 8 Slide the retaining clip (1) up to lock the blade adapter in place.

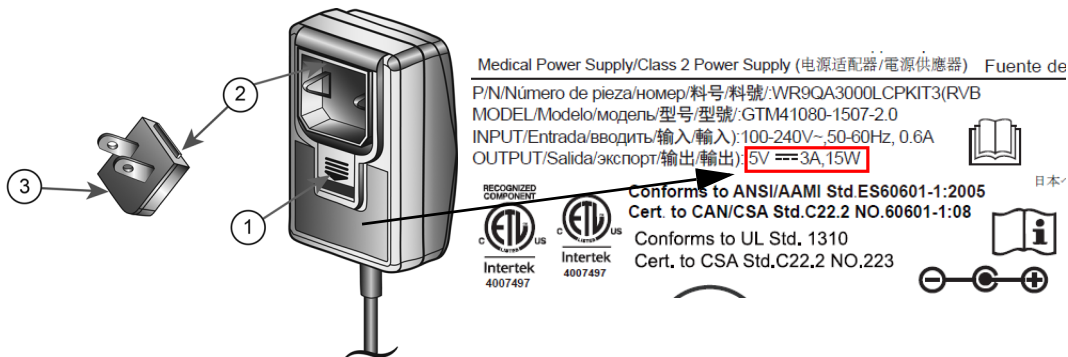CAUTION

THE ANALYZER WILL NOT CHARGE IF THE POWER ADAPTER OUTPUT IS GREATER THAN 5V 3A, 15W. IF THE BATTERY ICON IS FLASHING OR THE BATTERY CHARGE IS REQUIRED, MAKE SURE THE ANALYZER IS CONNECTED TO THE BD SUPPLIED POWER ADAPTER.

- 9 Insert the coaxial power connector (identified by the BD label near the connector) into the round port on the back of the Analyzer, and insert the power adapter plug into an electrical outlet.
- 10 Perform a verification test as described in section 5.1.1.

2.4 Printer Setup – Optional

Printer Unpacking and Setup

The BD Veritor™ Plus Analyzer can be configured to print test results to a connected printer. The printed language will match the currently configured display language. The BD Veritor™ Plus Analyzer supports only specific printer models. Contact your BD representative to obtain a list of compatible printer models and operating instructions. A USB 2.0 On-The-Go micro-A to B cable, or USB On-The-Go adapter, is required. The cable may be purchased from computer supply retailers or may be ordered from BD under catalog number 443907.

For unpacking and setup instructions, or information about ordering replacement paper, refer to the instructions provided with the printer.

The printer is connected to the BD Veritor™ Plus Analyzer as follows:

- 1 Connect the USB B end of the USB cable to the corresponding port on the back of the printer.
- 2 Connect the USB micro-A end of the USB cable to the corresponding rectangular port on the back of the BD Veritor™ Plus Analyzer.

NOTE

The BD Veritor™ Plus Analyzer will only print test results when the printer is connected and powered on during the test procedure. Due to data security considerations, test results cannot be recalled for subsequent printing. Prior to initiating testing, ensure that the printer paper supply is adequate to support the test volume.

2.5 Module and Accessory Specifications – Optional

Optional Modules and Spare Parts	BD Catalog Number
BD Veritor™ InfoScan	256068
BD Veritor™ InfoWiFi module	445010
BD Veritor™ Plus Connect	444881
BD Veritor™ Plus Analyzer Verification Cartridge	256047
BD Veritor™ Plus AC Power Adapter	445201
BD Veritor™ System USB Cable	445202
BD Veritor™ System Printer Cable	443907
Module Environmental Specifications	
Operating Conditions	
Ambient Temperature	15.0–30.0 °C (59–86 °F)
Ambient Humidity	Up to 85% RH, non-condensing
Ambient Light	0–3,000 lux, uniform illumination
Ambient Environment	Indoor Use Only

Other	
Overvoltage Category	II (per IEC 60664)
Lifetime Specification	The BD Veritor™ barcode reader equipped modules do not expire. The BD Veritor™ Plus Connect does not expire.
Performance Specifications	
Compatible 1D Barcode Protocol	Code 128, Code 39, UPC / EAN, Interleaved 2 of 5, Codabar, GS1 DataBar, GS1-128. ASCII-128 character set.
Operator ID Maximum Barcode Length	20 characters (Longer barcodes are truncated to 20 characters)
Specimen ID Maximum Barcode Length	20 characters (Longer barcodes are truncated to 20 characters)
Data Transmission Security Protocols	HTTPS with SSL/TLS certification. AES-256 payload encryption.

Default Workflow Configuration	
Operator Identification Scan	Disabled
Specimen Identification Scan	Enabled
Kit Lot Information Scan	Disabled
Display Language	English
USB Port	Locked

2.6 BD Veritor™ InfoScan Module Setup – Optional

Addition of the BD Veritor™ InfoScan module enables the Analyzer to capture barcoded information and integrate it into the test record. Utilized together with BD Veritor™ Plus Connect, and the BD Synapsys™ Informatics Solutions it allows secure transmission of test results to a facility LIS. The included USB cable may be used to network the Analyzer either to a printer or to the BD Veritor™ Plus Connect.

Steps for installing the module are:

- 1 Turn the Analyzer power off.
- 2 Turn the Analyzer upside down.

CAUTION

**Ensure that the power button is not accidentally pressed while the Analyzer is upside down.
Power to the Analyzer should remain off during module installation.**

- 3 Depress the locking tab on the module bay cover.

WARNING

ELECTRICAL EXPOSURE HAZARD. DO NOT TOUCH THE ELECTRICAL CONTACTS INSIDE THE MODULE BAY.

-
- 4 Slide the module bay cover off of the Analyzer. The module bay cover should be stored in a secure location.
 - 5 Remove protective plastic cover from the module interface pins if present.
 - 6 Align the module alignment features to the Analyzer track.
 - 7 Slide the module into the Analyzer until an audible click is heard.
 - 8 The Analyzer is now ready for operation. See section 4.6 for optional configuration instructions.

2.7 BD Veritor™ Plus Connect Set Up – Optional

The BD Veritor™ Plus Connect accessory utilized together with a barcode reader equipped module allows secure transmission of test results to a facility LIS. Additional IT setup is required in order to install and utilize the BD Veritor™ Plus Connect accessory. To create or troubleshoot this connection, please contact your local BD representative. Installation instructions for site IT staff are further detailed in documentation included with BD Veritor™ Plus Connect.

2.8 BD Veritor™ InfoWiFi Module Setup – Optional

The addition of the BD Veritor™ InfoWiFi module enables the Analyzer to capture barcoded information and integrate it into the test record. It then allows secure transmission of test results to a facility LIS. The USB cable included with the BD Veritor™ InfoWiFi module can be used to connect the Analyzer to the BD Veritor™ Plus Connect. In order to complete the connection to a wireless network, additional IT setup is required. Installation instructions for site IT staff are further detailed in documentation included with BD Veritor™ InfoWiFi.

To install InfoWiFi module:

- 1 Turn off the Analyzer power.
- 2 Turn the Analyzer upside down.

CAUTION

**Ensure that the power button is not accidentally pressed while the Analyzer is upside down.
Power to the Analyzer should remain off during module installation.**

- 3 Depress the locking tab on the module bay cover.

WARNING

ELECTRICAL EXPOSURE HAZARD. DO NOT TOUCH THE ELECTRICAL CONTACTS INSIDE THE MODULE BAY.

- 4 Slide the module bay cover off of the Analyzer. The module bay cover should be stored in a secure location.
- 5 Remove protective plastic cover from the module interface pins if present.
- 6 Align the module alignment features to the Analyzer track.
- 7 Slide the module into the Analyzer until an audible click is heard.
- 8 The Analyzer is now ready for operation. See section 4.6 for optional configuration instructions.

NOTE

The BD Veritor™ Plus InfoWiFi module does not expire and can be transferred to a different BD Veritor™ Plus Analyzer. Additional setup and registration steps may be required. See documentation provided with the BD Veritor™ InfoWiFi module.

- 9 See instructions included with the BD Veritor™ InfoWiFi module for configuration details and associated barcodes.

3 Controls and Indicators

3.1 General

This section describes the meaning and use of the controls and indicators of the BD Veritor™ Plus Analyzer. Some components are illustrated in the figures accompanying the related text.

The following controls and indicators are discussed:

- Coaxial Power Port
- Power Button
- System Indicators
- Test Device Slot
- USB Port

WARNING
**ALL USERS SHOULD BECOME THOROUGHLY FAMILIAR WITH ALL CONTROLS AND INDICATORS
BEFORE ATTEMPTING TO OPERATE THE INSTRUMENT.**

3.1.1 Coaxial Power Port

The BD Veritor™ Plus Analyzer features a rechargeable lithium-ion battery for portable usage and includes an AC power adapter to recharge the battery or to directly power the unit from a facility power supply.

To charge the BD Veritor™ Plus Analyzer, plug the cylindrical end of the charger into the corresponding port on the back of the Analyzer and plug the power adapter into the facility power source.

The BD Veritor™ Plus Analyzer can be externally powered only from the coaxial power port on the rear of the instrument. The Analyzer cannot be powered through its USB port. When the BD Veritor™ Plus Analyzer is properly connected to the facility power supply and powered on, the Plug icon will appear on the screen as shown below in section 3.1.3.

3.1.2 Power Button

To power on the Analyzer, depress the power button. After doing so, the BD Veritor™ Plus Analyzer screen will display the BD Logo, then enter a self-test mode. After the self-test is completed successfully, the Analyzer is ready to test and will prompt the user to insert a test device.

BD Veritor™ Plus Analyzers that have an optional module installed may prompt the user for additional workflow steps during the power-on process. For instructions relating to these steps, see section 4.

To power off the Analyzer, depress the power button for at least half a second and release.

If the Analyzer is left unattended for 15 minutes (when operating on the internal battery), or 60 minutes (when operating on the external power adapter), the Analyzer will automatically shut off and the test result will not be retained on the screen.


3.1.3 System Indicators

The meanings of on-screen displays used on the BD Veritor™ Plus Analyzer are as follows:

Symbol	Meaning
	<p>Hourglass: An indication that the Analyzer is nearing end of life, either due to usage duration (time) or usage count. In addition to the Hourglass icon, a descriptive message will appear on the screen after the power-on self-test is complete. The message describes which end-of-life event is approaching and the remaining usage time or usage count available. Each end-of-life event generates a unique message; insertion of a test device terminates all messages and immediately initiates a test.</p> <p>When this symbol is present, the user should contact their BD representative to order a new BD Veritor™ Plus Analyzer.</p> <p>The symbol switches from a continuous state to flashing when end of life is imminent.</p> <p>See section 2.2 for an explanation of lifetime specifications.</p>
	<p>Plug: When the Analyzer is properly connected to the facility power supply and powered on, the Plug icon will appear on screen.</p> <p>Before calling BD, check if the supplied AC Adapter is used.</p> <p>Battery Not Charging:</p> <p>In the event the battery icon does not indicate the battery is charging while the power supply is plugged into the BD Veritor™ Plus Analyzer, examine the AC Adapter plugged into the Analyzer. The supplied AC Adapter has a BD label near the coaxial power connector and conforms to the sample labeling in section 2.3.</p> <p>If the Plug icon is flashing, it indicates a failure in the battery charging circuit. Contact your local BD representative. The Analyzer may continue to be used when connected to the facility power supply.</p>
	<p>Battery: An indication of battery strength. The maximum battery strength is represented by four bars.</p> <p>When the battery is nearly discharged, no bars are displayed and the battery outline flashes, indicating an imminent need for recharge. While the battery is charging, the bars will cycle to indicate charging is underway. Once the battery is fully charged, the cycling will stop and four bars will appear in a continuous state.</p>
	<p>On BD Veritor™ Plus Analyzers that are connected to the BD Veritor™ Plus Connect, this symbol indicates when the Analyzer is connected to the facility LIS.</p> <p>On BD Veritor™ Plus Analyzers equipped with a BD Veritor™ InfoWiFi module, this symbol indicates connectivity to the wireless network. Please refer to the instructions provided with the module for more details.</p>
	<p>This symbol only appears on BD Veritor™ Plus Analyzers that are equipped with the BD Veritor™ connectivity options. The envelope symbol appears when test results are waiting to be transmitted. A flashing envelope indicates that data are in the process of being transmitted. The envelope symbol disappears when result transmission is complete.</p>
	<p>Do not use if product damaged - Indicates a medical device that should not be used if the package has been damaged or opened.</p>
	<p>Disposal of electronic products / Separate collection for waste of electrical and electronic equipment. Do not dispose of battery in municipal waste.</p>

3.1.4 Test Device Slot

The test device slot is located on the right side of the BD Veritor™ Plus Analyzer (see arrow in figure below). Test devices are inserted by aligning the insertion arrow on the device with the insertion arrow on the Analyzer. A distinctive “click” will be noted when the test device is properly engaged in the Analyzer.

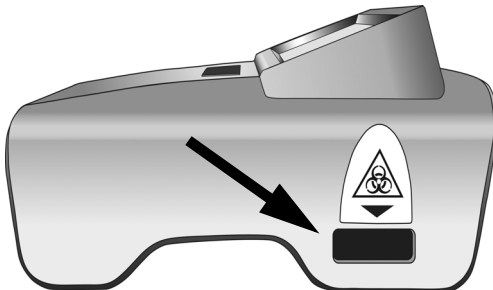3.1.5 USB Port

The BD Veritor™ Plus Analyzer includes a USB On-The-Go port which is used for four applications:

- Analyzer Update: the BD Veritor™ Plus Analyzer can be updated, using a USB flash drive that is provided by BD, to add capability to read additional BD Veritor™ assays and to add new features. Instructions to perform these updates are provided with the USB flash drive when such updates are available.
- Printing: An optional printer is connected to the BD Veritor™ Plus Analyzer via this port. See section 2.4 for instructions.
- Data Download: BD Veritor™ Plus Analyzers that are equipped with a barcode scanner can be unlocked to enable secure transfer of the test record to a computer via this port. See section 4.8 for instructions.
- Data Transmission: BD Veritor™ Plus Analyzers that are equipped with the BD Veritor™ Plus Connect accessory can transmit test results to a facility LIS over a secure connection via this port.

NOTE

The BD Veritor™ Plus Analyzer can be externally powered only from the coaxial power port on the rear of the instrument. The Analyzer cannot be powered through its USB port.

4 Operation

4.1 General

Section 4.2 provides instructions for routine operation of the BD Veritor™ Plus Analyzer using the Analyze Now operational mode. In this mode, the test device is prepared and developed outside of the Analyzer, then inserted to perform the final read. Section 4.3 provides instructions for routine operation of the BD Veritor™ Plus Analyzer using the Walk Away operational mode. In this mode, the test device is prepared outside the Analyzer but allowed to develop inside the Analyzer. The Analyzer initiates the final read at the appropriate time. Both sections include instructions for operation with or without a BD Veritor™ barcode reader equipped module installed. Sections 4.4 and 4.5 describe the Quality Control (QC) test procedure in Analyze Now and Walk Away mode, respectively.

4.2 Using the BD Veritor™ Plus Analyzer – Analyze Now Operation

The BD Veritor Plus Analyzer is intended for use with BD Veritor™ Test Devices. Specific instructions for performing a BD Veritor™ Test can be found in the package insert included with each test kit. After following those instructions the procedure for interpreting a test result with the BD Veritor™ Plus Analyzer using Analyze Now operation is as follows:

- 1 Place the BD Veritor™ Plus Analyzer on a flat, dry, stable surface.
- 2 Ensure that the BD Veritor™ Plus Analyzer is not in direct sunlight or exposed to a bright light.
- 3 Check the BD Veritor™ Plus Analyzer for an inserted device. If a device is present, remove it from the Analyzer.

-
- 4 Press the front panel power button. The BD Veritor™ Plus Analyzer will complete a self-test before it is ready for use. After the self-test completes and any temporary messages are presented, the display window shows INSERT TEST DEVICE OR DOUBLE-CLICK BUTTON FOR WALK AWAY MODE.

If a BD Veritor™ barcode reader equipped module is installed:

The display window shows SCAN CONFIG BARCODE for 2 seconds after the self-test. This is an optional step, no action is required. To perform instrument configuration, refer to section 4.6 for configuration instructions.

- 5 When the display window shows INSERT TEST DEVICE OR DOUBLE-CLICK BUTTON FOR WALK AWAY MODE, insert the fully developed BD Veritor™ Test Device into the slot on the right side of the Analyzer, aligning the insertion arrow on the test device with the arrow above the slot. Insert the device fully until it stops. A distinctive “click” will be noted when the device is properly aligned in the Analyzer.

If a BD Veritor™ barcode reader equipped module is installed:

- a If the BD Veritor™ Plus Analyzer has been configured with barcode ID – Enable Operator ID, it will display SCAN OPERATOR ID. When this message is displayed, after verifying that the Operator ID to be scanned is correct, the operator should scan his/her identification barcode. This message is displayed for up to 30 seconds. If scans are not completed in this 30 second period, the BD Veritor™ Test Device must be removed and the read process restarted with step 5. After a specific Operator Identification has been recorded, that Operator Identification will be utilized for all subsequent tests until the BD Veritor™ Plus Analyzer is powered off. On the next power cycle, the SCAN OPERATOR ID prompt will be repeated during the first test. A test result cannot be generated if an Operator ID is not scanned when the Operator ID workflow option is enabled.

BD recommends reading Operator Identification barcodes with the Analyzer positioned at the edge of a flat surface. The barcode should then be moved toward the barcode window to be scanned. The scanned barcode value will be shown in the next display window.

- b If the BD Veritor™ Plus Analyzer has been configured with barcode ID – Enable Specimen ID, it will display SCAN SPECIMEN ID. When this message is displayed, after verifying that the Specimen ID is correct, the operator should scan the specimen's identification barcode. This message is displayed for up to 30 seconds. If scans are not completed in this 30 second period, the BD Veritor Test Device must be removed and the read process restarted with step 5. The Specimen ID scan prompt is repeated for every test. A test result cannot be generated if a Specimen ID is not scanned when the Specimen ID workflow option is enabled.

BD recommends reading Specimen Identification barcodes with the Analyzer positioned at the edge of a flat surface. The barcode should then be moved towards the barcode window to be scanned. The scanned barcode value will be shown in the next display window.

- c If the BD Veritor™ Plus Analyzer has been configured with barcode ID – Enable Kit Lot Information, it will display SCAN KIT LOT NUMBER. When this message is displayed, the operator should scan the barcode on the exterior of the BD Veritor™ Test Kit box. This message is displayed for up to 30 seconds. If scans are not completed in this 30 second period, the BD Veritor™ Test Device must be removed and the read process restarted with step 5. The Kit Lot Number scan prompt is repeated for every test. A test result cannot be generated if the Kit Lot Number is not scanned when the Kit Information workflow option is enabled.

If the BD Veritor™ Test Kit box label has two barcodes, scan only the upper barcode which starts with (17).

BD recommends reading the Kit Information barcode with the Analyzer positioned at the edge of a flat surface. The barcode should then be moved toward the barcode window to be scanned. The scanned barcode value will be shown in the next display window.

- 6 After insertion of the test device, the BD Veritor™ Plus Analyzer will progress through two processing steps: a reading step, followed by an analyzing step. The display will show the remaining time for each step as they are performed. Do not touch the Analyzer or remove the test device during this time.
- 7 When the analysis is complete, the test result will be displayed with the name of the test and a result. Refer to the package insert for the BD Veritor™ Test Device for specific result interpretations. If a printer is connected to the unit and powered on, the test result will automatically be sent to the printer.

If a BD Veritor™ barcode reader equipped module is installed:

If the BD Veritor™ Plus Analyzer has been configured with barcode ID – Enable Specimen ID, the specimen identification is also displayed on the screen.

- 8 Verify that the test type and specimen identification (if appropriate) are correct.
- 9 Once the test type and specimen identification, if appropriate, is verified and the result noted, remove the test device by pulling it out. The display will show INSERT TEST DEVICE OR DOUBLE-CLICK BUTTON FOR WALK AWAY MODE to indicate the Analyzer is ready to perform another test.

If a BD Veritor™ Plus Connect accessory or a connectivity module is installed:

The ENVELOPE symbol will appear to indicate that results are being transmitted. In the event that the BD Veritor™ Plus Analyzer cannot transmit the results to BD Synapsys™ Informatics Solution, it will queue all results to be transmitted and continuously attempt to transmit the results while it is powered on. If the Analyzer is powered off while the ENVELOPE symbol is still present, the Analyzer will queue the result and transmit it the next time that it is powered up. The symbol will disappear after the results have been transmitted.

- 10 To initiate a new test, repeat steps 5–9 above, or turn off the power by depressing the front panel power button for at least half a second and releasing it. If the Analyzer is left unattended for 15 minutes (when operating on the internal battery), or 60 minutes (when operating on the external power adapter), the BD Veritor™ Plus Analyzer will automatically shut off and the test result will not be retained on the screen.

4.3 Using the BD Veritor Plus Analyzer – Walk Away Operation

NOTE

The AC power adapter must be connected to the Analyzer and plugged into a facility power source to use Walk Away operation.

NOTE

If the Analyzer is left unattended for 60 minutes after completion of the test run, the power will automatically shut off and the test result will not be retained on the screen.

The BD Veritor™ Plus Analyzer is intended for use with BD Veritor™ Test Devices. Specific instructions for performing a BD Veritor™ Test can be found in the package insert included with each test kit. After following those instructions up to sample preparation in the reagent tube, the procedure for interpreting a test result with the BD Veritor™ Plus Analyzer using Walk Away operation is as follows. DO NOT DISPENSE THE SPECIMEN INTO THE TEST DEVICE SAMPLE WELL AT THIS TIME.

- 1 Place the BD Veritor™ Plus Analyzer on a flat, dry, stable surface.
- 2 Ensure that the BD Veritor™ Plus Analyzer is not in direct sunlight or exposed to a bright light.
- 3 Ensure that the AC power adapter is connected to the Analyzer and plugged into a facility power source.
- 4 Check the BD Veritor™ Plus Analyzer for an inserted device. If a device is present, remove it from the Analyzer.
- 5 Press the front panel power button. The BD Veritor™ Plus Analyzer will complete a self-test before it is ready for use. After the self-test completes and any temporary messages are presented, the display window shows INSERT TEST DEVICE OR DOUBLE-CLICK BUTTON FOR WALK AWAY MODE.

If a BD Veritor™ barcode reader equipped module is installed:

The display window shows SCAN CONFIG BARCODE for 2 seconds after the self-test. This is an optional step, no action is required. To perform instrument configuration, refer to section 4.6 for configuration instructions.

- 6 When the display window shows INSERT TEST DEVICE OR DOUBLE-CLICK BUTTON FOR WALK AWAY MODE, double-click the power button.

If a BD Veritor™ barcode reader equipped module is installed:

- a If the BD Veritor™ Plus Analyzer has been configured with barcode ID – Enable Operator ID, it will display SCAN OPERATOR ID. When this message is displayed, after verifying that the Operator ID to be scanned is correct, the operator should scan his/her identification barcode. This message is displayed for up to 30 seconds. If scans are not completed in this 30 second period, the process will restart with step 6. After a specific Operator Identification has been recorded, that Operator Identification will be utilized for all subsequent tests until the BD Veritor™ Plus Analyzer is powered off. On the next power cycle, the SCAN OPERATOR ID prompt will be repeated during the first test. A test result cannot be generated if the Operator ID is not scanned when the Operator ID workflow option is enabled.

BD recommends reading Operator Identification barcodes with the Analyzer positioned at the edge of a flat surface. The barcode should then be moved toward the barcode window to be scanned. The scanned barcode value will be shown in the next display window.

-
- b If a BD Veritor™ Plus Analyzer has been configured with barcode ID – Enable Specimen ID, it will display SCAN SPECIMEN ID. When this message is displayed, after verifying that the Specimen ID is correct, the operator should scan the specimen's identification barcode. This message is displayed for up to 30 seconds. If scans are not completed in this 30 second period, the process will restart with step 6. The Specimen ID scan prompt is repeated for every test. A test result cannot be generated if the Specimen ID is not scanned when the Specimen ID workflow option is enabled.

BD recommends reading Specimen Identification barcodes with the Analyzer positioned at the edge of a flat surface. The barcode should then be moved toward the barcode window to be scanned. The scanned barcode value will be shown in the next display window.
 - c If the BD Veritor™ Plus Analyzer has been configured with barcode ID – Enable Kit Lot Information, it will display SCAN KIT LOT NUMBER. When this message is displayed, the operator should scan the barcode on the exterior of the BD Veritor™ Test Kit box. This message is displayed for up to 30 seconds. If scans are not completed in this 30 second period, the process restarts with step 6. The Kit Lot Number scan prompt is repeated for every test. A test result cannot be generated if the Kit Lot Number is not scanned when the Kit Information workflow option is enabled.

If the BD Veritor™ Test Kit box label has two barcodes, scan only the upper barcode which starts with (17).

BD recommends reading the Kit Information barcode with the Analyzer positioned at the edge of a flat surface. The barcode should then be moved toward the barcode window to be scanned. The scanned barcode value will be shown in the next display window.
 - 7 The display window will now show ADD SPECIMEN TO TEST DEVICE AND INSERT IMMEDIATELY. This message is displayed for up to 3 minutes, after which the process restarts with step 6. Apply the prepared specimen to the test device sample well as instructed in the test device kit package insert and immediately insert the test device into the slot on the right side of the Analyzer, aligning the insertion arrow on the test device with the arrow above the slot. Insert the device fully until it stops. A distinctive "click" will be noted when the device is properly aligned in the Analyzer. During this process the test device must remain horizontal to prevent spilling the specimen out of the sample well.
 - 8 The display window will now show DO NOT DISTURB TEST IN PROGRESS. The incubation time is determined based on the test device barcode. The incubation time remaining is shown on the display. Do not disturb the test device or Analyzer during the incubation period. Do not remove the test device; doing so will cause the test to abort.
 - 9 After the incubation is complete, the BD Veritor™ Plus Analyzer will progress through two processing steps: a reading step, followed by an analyzing step. The display will show the remaining time for each step as they are performed. Do not touch the Analyzer or remove the test device during this time.
 - 10 When the analysis is complete, the test result will be displayed with the name of the test and a result. Refer to the package insert for the BD Veritor™ Test Device for specific result interpretations. If a printer is connected to the unit and powered on, the test result will automatically be sent to the printer.

If the BD Veritor™ barcode reader equipped module is installed:

If the BD Veritor™ Plus Analyzer has been configured with barcode ID – Enable Specimen ID, the specimen identification is also displayed on-screen.

- 11 Verify that the test type and specimen identification (if appropriate) are correct.
- 12 Once the test type and specimen identification, if appropriate, is verified and the result noted, remove the test device by pulling it out. The display will show INSERT TEST DEVICE OR DOUBLE-CLICK BUTTON FOR WALK AWAY MODE to indicate the Analyzer is ready to perform another test.

If a BD Veritor™ Plus Connect accessory or a connectivity module is installed:

The ENVELOPE symbol will appear to indicate that results are being transmitted. In the event that the BD Veritor™ Plus Analyzer cannot transmit the results to BD Synapsys™ Informatics Solution, it will queue all results to be transmitted and continuously attempt to transmit the results while it is powered on. If the Analyzer is powered off while the ENVELOPE symbol is still present, the Analyzer will queue the result and transmit it the next time that it is powered up. The symbol will disappear after the results have been transmitted.

- 13 To initiate a new test, repeat steps 6–12 above, or turn off the power by depressing the front panel power button for at least half a second and releasing it. If the Analyzer is left unattended for 60 minutes, the BD Veritor™ Plus Analyzer will automatically shut off and the test result will not be retained on the screen.

4.4 Kit Quality Control (QC) Test Procedure – Analyze Now Operation

NOTE

To utilize the Analyzer's QC Test documentation feature, the Analyzer must have a barcode reading module installed, and must have been configured with barcode – Enable Specimen ID.

NOTE

The Analyzer records the Kit Lot and Control Swab expiration dates in the test record but does not restrict the use of expired reagents. Management of expired materials is the responsibility of the user per local site regulations or practices. BD recommends that expired materials never be used.

The BD Veritor™ Plus Analyzer is intended for use with BD Veritor™ Test Devices. Specific instructions for performing a BD Veritor™ Quality Control (QC) Test can be found in the package insert included with each test kit. After following those instructions the procedure for interpreting a test result with the BD Veritor™ Plus Analyzer using Analyze Now operation is as follows:

- 1 Place the BD Veritor™ Plus Analyzer on a flat, dry, stable surface.
- 2 Ensure that the BD Veritor™ Plus Analyzer is not in direct sunlight or exposed to a bright light.
- 3 Check the BD Veritor™ Plus Analyzer for an inserted device. If a device is present, remove it from the Analyzer.
- 4 Press the front panel power button. The BD Veritor™ Plus Analyzer will complete a self-test before it is ready for use. After the self-test completes, the display window shows SCAN CONFIG BARCODE. This is an optional step, no action is required. After 2 seconds, the display window shows INSERT TEST DEVICE OR DOUBLE-CLICK BUTTON FOR WALK AWAY MODE. To perform instrument configuration, refer to section 4.6 for configuration instructions.
- 5 When the display window shows INSERT TEST DEVICE OR DOUBLE-CLICK BUTTON FOR WALK AWAY MODE, insert the fully developed BD Veritor™ Test Device which has been prepared with an external control into the slot on the right side of the Analyzer, aligning the insertion arrow on the test device with the arrow above the slot. Insert the device fully until it stops. A distinctive "click" will be noted when the device is properly aligned in the Analyzer.
- 6 If the BD Veritor™ Plus Analyzer has been configured with barcode ID – Enable Operator ID, it will display SCAN OPERATOR ID. When this message is displayed, after verifying that the Operator ID to be scanned is correct, the operator should scan his/her identification barcode. This message is displayed for up to 30 seconds. If scans are not completed in this 30 second period, the BD Veritor™ Test Device must be removed and the read process restarted with step 5. After a specific Operator Identification has been recorded, that Operator Identification will be utilized for all subsequent tests until the BD Veritor™ Plus Analyzer is powered off. On the next power cycle, the SCAN OPERATOR ID prompt will be repeated during the first test. A test result cannot be generated if the Operator ID is not scanned when the Operator ID workflow option is enabled.

BD recommends reading Operator Identification barcodes with the BD Veritor™ Plus Analyzer positioned at the edge of a flat surface. The barcode should then be moved toward the barcode window to be scanned. The scanned barcode value will be shown in the next display window.

- 7 **In the case of Quality Control (QC) Tests, the barcode on the printed label of the swab control supplied with the kit is used as the Specimen ID.**

The display window will now display SCAN SPECIMEN ID. When this message is displayed, the operator should scan the barcode on the printed label on the swab control supplied with the kit. This message is displayed for up to 30 seconds. If scans are not completed in this 30 second period, the BD Veritor™ Test Device must be removed and the read process restarted with step 5. A test result cannot be generated if the swab control barcode is not scanned.

BD recommends reading Quality Control (QC) barcodes with the Analyzer positioned at the edge of a flat surface. The barcode should then be moved toward the barcode window to be scanned. The scanned barcode value will be shown in the next display window.

- 8 **In the case of Quality Control (QC) Tests, the barcode on the printed label of the test device kit box must be scanned.**

The display window will show SCAN KIT LOT NUMBER whether the kit information workflow configuration setting is enabled or disabled for normal operation. When this message is displayed, the operator should scan the barcode on the exterior of the BD Veritor™ Test Kit box. This message is displayed for up to 30 seconds. If scans are not completed in this 30 second period, the Test Device must be removed and the read process restarted with step 5. A test result cannot be generated if the Kit Lot Number is not scanned.

If the BD Veritor™ Test Kit box label has two barcodes, scan only the upper barcode which starts with (17).

BD recommends reading Kit Lot Information barcodes with the BD Veritor™ Plus Analyzer positioned at the edge of a flat surface. The barcode should then be moved toward the barcode window to be scanned. The scanned barcode value will be shown in the next display window.

- 9 After all required scans are complete, the BD Veritor™ Plus Analyzer will progress through two processing steps: a reading step, followed by an analyzing step. The display will show the remaining time for each step as they are performed. Do not touch the Analyzer or remove the test device during this time.

- 10 When the analysis is complete, the test result will be displayed with the control swab identification, the name of the test, a result, and an overall "QC PASS" or "QC FAIL" indication. Refer to the package insert for the BD Veritor™ Test Device for specific result interpretations. If a printer is connected to the unit and powered on, the test result will automatically be sent to the printer.
- 11 Verify that the test type and QC swab control are correct.
- 12 Once the test type and swab control are verified and the result noted, remove the test device by pulling it out. The display will show INSERT TEST DEVICE OR DOUBLE-CLICK BUTTON FOR WALK AWAY MODE to indicate the Analyzer is ready to perform another test.
- 13 When the BD Veritor™ Plus Connect accessory is installed, the ENVELOPE symbol will appear to indicate that results are being transmitted over the secure connection. In the event that the BD Veritor™ Plus Analyzer cannot transmit the results to BD Synapsys™ Informatics Solution, it will queue all results to be transmitted and continuously attempt to transmit the results while it is powered on. If the Analyzer is powered off while the ENVELOPE symbol is still present, the Analyzer will queue the result and transmit it the next time that it is powered up. The symbol will disappear after the results have been transmitted.
- 14 To initiate a new test, repeat steps 5–13 above, or turn off the power by depressing the front panel power button for at least half a second and releasing it. If the Analyzer is left unattended for 15 minutes (when operating on the internal battery), or 60 minutes (when operating on the external power adapter), the BD Veritor™ Plus Analyzer will automatically shut off and the test result will not be retained on the screen.

4.5 Kit Quality Control (QC) Test Procedure – Walk Away Operation

NOTE

To utilize the Analyzer's QC Test documentation feature, the Analyzer must have a barcode reading module installed and must have been configured with barcode–Enable Specimen ID.

NOTE

The Analyzer records the Kit Lot and Control Swab expiration dates in the test record but does not restrict the use of expired reagents. Management of expired materials is the responsibility of the user per local site regulations or practices. BD recommends that expired materials never be used.

NOTE

The AC power adapter must be connected to the Analyzer and plugged into a facility power source to use Walk Away operation.

NOTE

If the Analyzer is left unattended for 60 minutes after completion of the test run, the power will automatically shut off and the test result will not be retained on the screen.

The BD Veritor™ Plus Analyzer is intended for use with BD Veritor™ Test Devices. Specific instructions for performing a BD Veritor™ Quality Control (QC) Test can be found in the package insert included with each test kit. After following those instructions up to sample preparation in the reagent tube, the procedure for interpreting a test result with the BD Veritor™ Plus Analyzer using Walk Away operation is as follows. DO NOT DISPENSE THE CONTROL INTO THE TEST DEVICE SAMPLE WELL AT THIS TIME.

- 1 Place the BD Veritor™ Plus Analyzer on a flat, dry, stable surface.
- 2 Ensure that the BD Veritor™ Plus Analyzer is not in direct sunlight or exposed to a bright light.
- 3 Ensure that the AC power adapter is connected to the Analyzer and plugged into a facility power source.
- 4 Check the BD Veritor™ Plus Analyzer for an inserted device. If a device is present, remove it from the Analyzer.
- 5 Press the front panel power button. The BD Veritor™ Plus Analyzer will complete a self-test before it is ready for use. After the self-test completes, the display window shows SCAN CONFIG BARCODE. This is an optional step, no action is required. After 2 seconds, the display window shows INSERT TEST DEVICE OR DOUBLE-CLICK BUTTON FOR WALK AWAY MODE. To perform instrument configuration, refer to section 4.6 for configuration instructions.

- 6 When the display window shows INSERT TEST DEVICE OR DOUBLE-CLICK BUTTON FOR WALK AWAY MODE, double-click the power button.
- 7 If the BD Veritor™ Plus Analyzer has been configured with barcode ID – Enable Operator ID, it will display SCAN OPERATOR ID. When this message is displayed, after verifying that the Operator ID to be scanned is correct, the operator should scan his/her identification barcode. This message is displayed for up to 30 seconds. If scans are not completed in this 30 second period, the process will restart with step 6. After a specific Operator Identification has been recorded, that Operator Identification will be utilized for all subsequent tests until the BD Veritor™ Plus Analyzer is powered off. On the next power cycle, the SCAN OPERATOR ID prompt will be repeated during the first test. A test result cannot be generated if the Operator ID is not scanned when the Operator ID workflow option is enabled.

BD recommends reading Operator Identification barcodes with the BD Veritor™ Plus Analyzer positioned at the edge of a flat surface. The barcode should then be moved toward the barcode window to be scanned. The scanned barcode value will be shown in the next display window.
- 8 **In the case of Quality Control (QC) Tests, the barcode on the printed label of the swab control supplied with the kit is used as the Specimen ID.**

The display window will now display SCAN SPECIMEN ID. When this message is displayed, the operator should scan the barcode on the printed label on the swab control supplied with the kit. This message is displayed for up to 30 seconds. If scans are not completed in this 30 second period, the process restarts with step 6. A test result cannot be generated if the swab control barcode is not scanned.

BD recommends reading Quality Control (QC) barcodes with the BD Veritor™ Plus Analyzer positioned at the edge of a flat surface. The barcode should then be moved toward the barcode window to be scanned. The scanned barcode value will be shown in the next display window.
- 9 **In the case of Quality Control (QC) Tests, the barcode on the printed label of the test device kit box must be scanned.**

The display window will show SCAN KIT LOT NUMBER whether the kit information workflow configuration setting is enabled or disabled for normal operation. When this message is displayed, the operator should scan the barcode on the exterior of the BD Veritor™ Test Kit box. This message is displayed for up to 30 seconds. If scans are not completed in this 30 second period, the process restarts with step 6. A test result cannot be generated if the Kit Lot Number is not scanned.

If the BD Veritor™ Test Kit box label has two barcodes, scan only the upper barcode which starts with (17).

BD recommends reading Kit Lot Information barcodes with the BD Veritor™ Plus Analyzer positioned at the edge of a flat surface. The barcode should then be moved toward the barcode window to be scanned. The scanned barcode value will be shown in the next display window.
- 10 The display window will now show ADD CONTROL TO TEST DEVICE AND INSERT IMMEDIATELY. This message is displayed for up to 3 minutes, after which the process restarts with step 6. Apply the prepared control to the test device sample well as instructed in the test device kit package insert and immediately insert the test device into the slot on the right side of the Analyzer, aligning the insertion arrow on the test device with the arrow above the slot. Insert the device fully until it stops. A distinctive “click” will be noted when the device is properly aligned in the Analyzer. During this process the test device must remain horizontal to prevent spilling the specimen out of the sample well.
- 11 The display window will now show DO NOT DISTURB TEST IN PROGRESS. The incubation time is determined based on the test device barcode. The incubation time remaining is shown on the display. Do not disturb the test device or Analyzer during the incubation period. Do not remove the test device; doing so will cause the test to abort.
- 12 After the incubation is complete, the BD Veritor™ Plus Analyzer will progress through two processing steps: a reading step, followed by an analyzing step. The display will show the remaining time for each step as they are performed. Do not touch the Analyzer or remove the test device during this time.
- 13 When the analysis is complete, the test result will be displayed with the control swab identification, the name of the test, a result, and an overall “QC PASS” or “QC FAIL” indication. Refer to the package insert for the BD Veritor™ Test Device for specific result interpretations. If a printer is connected to the unit and powered on, the test result will automatically be sent to the printer.
- 14 Verify that the test type and QC swab control are correct.
- 15 Once the test type and swab control are verified and the result noted, remove the test device by pulling it out. The display will show INSERT TEST DEVICE OR DOUBLE-CLICK BUTTON FOR WALK AWAY MODE to indicate the Analyzer is ready to perform another test.
- 16 When the BD Veritor™ Plus Connect accessory is installed, the ENVELOPE symbol will appear to indicate that results are being transmitted over the secure connection. In the event that the BD Veritor™ Plus Analyzer cannot transmit the results to BD Synapsys™ Informatics Solution, it will queue all results to be transmitted and continuously attempt to transmit the results while it is powered on. If the Analyzer is powered off while the ENVELOPE symbol is still present, the Analyzer will queue the result and transmit it the next time that it is powered up. The symbol will disappear after the results have been transmitted.
- 17 To initiate a new test, repeat steps 6–16 above, or turn off the power by depressing the front panel power button for at least half a second and releasing it. If the Analyzer is left unattended for 60 minutes, the BD Veritor™ Plus Analyzer will automatically shut off and the test result will not be retained on the screen.

4.6 Workflow Configuration Procedure – Optional

This section describes how to use the features of the BD Veritor™ Plus Analyzer barcode scanning module to configure various workflow options, such as selecting whether to require scanning an Operator ID or Specimen ID each time a test is performed. The scanner may also be used to enter various maintenance modes. See section 6 for more information on maintenance modes.

The Analyzer default settings appear in section 2.5.

To set workflow options, the barcode scanner turns on for 2 seconds every time the Analyzer is powered on. If a configuration barcode is not detected within those 2 seconds, the normal power up sequence continues. As the configuration scan duration is only 2 seconds, the user should have the desired configuration barcode preselected from pages 38–39 and ready to present to the Analyzer before turning on the Analyzer power.

The workflow configuration procedure is as follows:

- 1 Check the BD Veritor™ Plus Analyzer for an inserted device. If a device is present, remove it from the Analyzer.
- 2 Locate the appropriate barcode in section 7 to set the desired option.
- 3 Hold the Analyzer with the scanning window aimed near the appropriate barcode.
- 4 Press the front panel power button to turn on the power. During the self-test sequence, the scanner's red LED will flash briefly. Immediately after the self-test completes, the display window shows SCAN CONFIG BARCODE and the scanner's red LED turns on for 2 seconds.
- 5 Move the Analyzer closer or farther from the page until the red light is covering the desired barcode. The Analyzer will beep and display a confirmatory message when the barcode is read.
- 6 Additional settings may be made by scanning another configuration barcode within the next 2 seconds. As many settings as desired may be configured in this way. If the Analyzer times out and additional configurations are desired, simply turn the Analyzer power off by holding down the power button for at least half a second and releasing it. Power up again and repeat the workflow configuration procedure.
- 7 Except for USB Port Unlock, workflow options need only be set once. The Analyzer will retain the setting for its lifetime. Settings may be changed at any time by completing the workflow configuration procedure again.
- 8 The USB Port Unlock barcode is scanned during the configuration scan period to enable download of test results to a computer. The USB port remains unlocked for only the current power cycle. If the power is turned off, the USB Port Unlock barcode must be scanned again. Each Analyzer has a unique USB Unlock barcode supplied on an adhesive label with the Analyzer. The label may be affixed to this manual for safe keeping.

4.7 Workflow Configuration Barcodes

Scan the barcodes in section 7 to set the following workflow options.

On-Screen Language

The screen display language may be set to one of six choices. A confirmatory screen will be displayed in the current language. All subsequent display screens will be in the new language. If a printer is connected, the printed language will be the same as the on-screen language.

Operator Identification Scan

Set the Analyzer to prompt for Operator ID once per power cycle.

Specimen Identification Scan

Set the Analyzer to prompt for Specimen ID once per test.

Kit Lot Information Scan

Set the Analyzer to prompt for Kit Lot Information once per test.

NOTE

The Analyzer records the Kit Lot expiration date in the test record but does not restrict the use of expired reagents. Management of expired materials is the responsibility of the user per local site regulations or practices. BD recommends that expired materials never be used.

USB Port Unlock

Scan the USB Port Unlock barcode supplied with the Analyzer to enable download of test results to a computer. The port only remains unlocked for the current power cycle.

4.8 Data Download to a PC

The BD Veritor™ Plus Analyzer stores a data record of every test performed on a secured drive. Accessing the data requires:

- Installation of an optional BD Veritor™ barcode reader equipped module (see section 2.6 and 2.7)
- Access to the USB Port Unlock Barcode, included with the Analyzer
- A USB cable (USB-A to micro-B) (included with the BD Veritor™ barcode reader equipped Modules)
- A PC with an available USB 2.0 port (not included)

After installation of a BD Veritor™ barcode reader equipped module, results can be accessed by:

- 1 Place the BD Veritor™ Plus Analyzer on a flat, dry, stable surface.
- 2 Ensure that the BD Veritor™ Plus Analyzer is not in direct sunlight or exposed to a bright light.
- 3 Check the BD Veritor™ Plus Analyzer for an inserted device. If a device is present, remove it from the Analyzer.
- 4 Press the front panel power button. The BD Veritor™ Plus Analyzer will complete a self-test before it is ready for use.
- 5 When the display window reads SCAN CONFIG BARCODE, scan the USB Unlock Barcode. If the drive has been successfully unlocked, the display window will read USB PORT UNLOCKED.
- 6 Connect the BD Veritor™ Plus Analyzer to a PC by plugging the micro-B connector into the USB port of the BD Veritor™ Plus Analyzer and the A plug into a free USB port on the PC. The on-board memory will be associated as "BD Analyzer" with two unused drive letters on the attached PC; only one drive is accessible. Using a file utility such as Microsoft® Windows® Explorer, copy the file RECORD.CSV to the PC. This file may be opened in a spreadsheet program such as Microsoft® Excel or in a text editor. Note that the on-board drive is read-only and cannot be modified. Also note that all records are reported in Coordinated Universal Time (UTC), also known as Greenwich Mean Time (GMT).

NOTE

After the secured drive is unlocked, it will remain unlocked until the power is turned off. To access the drive again, the unlock barcode scan procedure must be repeated. Users should take care to cycle power immediately after completing the data download process to prevent unauthorized access.

5 Maintenance

5.1 General Procedures for Maintenance and Servicing

The BD Veritor™ Plus Analyzer requires little maintenance from the user to provide reliable performance. Any maintenance or repair not described in this section should be performed by BD representatives only.

WARNING

ALL MAINTENANCE AND REPAIR OTHER THAN THE PROCEDURES DESCRIBED IN THIS SECTION MUST BE PERFORMED BY QUALIFIED SERVICE PERSONNEL. NON-COMPLIANCE WITH THIS WARNING MAY RESULT IN PERSONAL INJURY OR INSTRUMENT MALFUNCTION.

5.1.1 BD Veritor™ System Verification Cartridge

NOTE

The BD Veritor™ System Verification Cartridge must be run using Analyze Now operation. Running in Walk Away mode is not allowed.

A Verification Cartridge is supplied to allow the user to perform a functional test on the Analyzer. The frequency of verification testing is dictated by local site regulations or practices. A verification test counts as one test towards the Analyzer's maximum allowed number of tests. All BD Veritor™ System Verification Cartridges must be obtained from BD or from a BD-authorized distributor. Cartridges from other manufacturers are not compatible with the BD Veritor™ Plus Analyzer.

To perform a verification test:

- 1 Place the BD Veritor™ Plus Analyzer on a flat, dry, stable surface.
- 2 Ensure that the BD Veritor™ Plus Analyzer is not in direct sunlight or exposed to a bright light.
- 3 Check the BD Veritor™ Plus Analyzer for an inserted device. If a device is present, remove it from the Analyzer.
- 4 Press the front panel power button. The BD Veritor™ Plus Analyzer must complete a self-test before it is ready for use.
- 5 When the display window reads INSERT TEST DEVICE OR DOUBLE-CLICK BUTTON FOR WALK AWAY MODE, remove the BD Veritor™ System Verification Cartridge from the pouch and insert into the slot on the right side of the Analyzer, aligning the insertion arrow on the Verification Cartridge with the arrow above the slot. Insert the cartridge fully until it stops. A distinctive "click" will be noted when the cartridge is properly engaged in the Analyzer.
- 6 The display will show the remaining read time as the read is performed. The read duration of the BD Veritor™ Plus Analyzer is approximately 30 seconds. Do not remove or touch the cartridge during this time.
- 7 When the read is complete, the test result will be displayed as one of the two following messages:
VERIFY PASS
or
VERIFY FAIL
CONSULT MANUAL
ERROR 25
- 8 If the message VERIFY PASS appears, the Analyzer may be used to test specimens. If the message VERIFY FAIL appears, the Analyzer may not be used to test specimens and may need to be replaced. See section 6.3 for troubleshooting information. Please contact your local BD representative for assistance.
- 9 Return the BD Veritor™ System Verification Cartridge to the resealable pouch for the next use.
- 10 If the verification cartridge is lost or damaged a replacement BD Veritor™ Plus Verification cartridge is available by ordering BD Catalog Number 256047. Please contact your local BD Representative for more information.

5.1.2 Cleaning

The outer case and display may be wiped with a clean towel lightly moistened with 70% isopropyl alcohol (IPA) or a 10% bleach solution. Do not introduce the cleaning solution or any other liquids directly into the unit. Do not use a saturated towel which may introduce liquid into the case or display seams. Ensure that the BD Veritor™ Plus Analyzer is dry and the surface is free of any residual cleaning solution prior to returning to use.

If a BD Veritor™ barcode reader equipped module is installed, BD does not recommend cleaning the barcode scanner window with any cleaning agent. Use a clean, soft towel lightly moistened with plain water to clean the window gently. Scratching the window may reduce the scanner's performance.

5.2 Servicing

There are no user serviceable components in the unit. For technical issues or questions, please contact your local BD representative.

5.3 Transportation and Storage

The Analyzer and optional modules should be transported in the original packaging at -20 °C – 65 °C for up to 2 days, -20 °C – 45 °C for up to 14 days, up to 85% RH, non-condensing. BD recommends the original packing materials be retained for this purpose. Follow local transportation regulations for shipping equipment containing a Li-ion secondary battery.

The BD Veritor™ Plus Analyzer should be stored at 15–30 °C, up to 85% RH, non-condensing.

6 Troubleshooting

6.1 Instrument Service

If your BD Veritor™ Plus Analyzer malfunctions or operates unusually in any way, you may initially attempt to solve the problem by following the recommendations in this section. All other servicing attempts will terminate the responsibility of the manufacturer under the terms of the warranty.

If you cannot repair a system malfunction, contact your local BD representative (contact numbers are listed on page 213).

6.2 Forced Power Down

If the BD Veritor™ Plus Analyzer is not responsive and does not turn off using the standard power down procedure, hold the power button down for approximately 10 seconds. The Analyzer will power down. Normal operation may then be resumed by following the standard power on sequence.


6.3 Troubleshooting Aids

The barcode scanner in the BD Veritor™ barcode reader equipped module may be used to initiate the display of system information or to enter a barcode scanner test mode.

Analyzer Information Display Mode

This mode will display the following system information:

- Analyzer serial number
- Analyzer expiration date and remaining test count
- Current configuration settings
- Firmware version

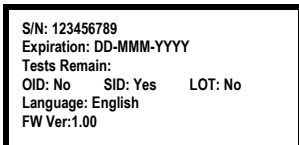To enter this mode, scan the following barcode during the workflow configuration scan step. The display will remain until the Analyzer power is turned off. See section 4.6 for information on the configuration scan step.

Analyzer Information Display

Show system information in the display window.

Barcode Scanner Test Mode

This mode places the BD Veritor™ barcode scanner into continuous operating mode. Any barcode presented to the scanner will be read and the value shown in the display window. This mode may be used to verify that the scanner is functioning correctly by scanning a known barcode, or may be used to confirm that the scanner is able to read the barcodes in use at a particular facility.

After scanning a barcode, its value will be shown in list form in the display window. Each value has an index number and the display window will scroll up to show new values. Barcode values longer than a single display line will wrap to the next line. If no barcode is presented within 30 seconds, the barcode scanner will turn off and a timeout message will appear in the display window. To exit the mode, turn off the power by depressing the front panel power button for at least half a second and releasing it. If the Analyzer is left unattended for 15 minutes (when operating on the internal battery) or 60 minutes (when operating on the external power adapter), the BD Veritor™ Plus Analyzer will automatically shut off and the display content will not be retained.

To enter this mode, scan the following barcode during the workflow configuration scan step. The test mode will remain active until the Analyzer power is turned off. See section 4.6 for information on the configuration scan step.

Barcode Scanner Test

Show scanned barcode values in the display window.

Battery Not Charging

In the event the battery icon does not indicate the battery is charging while the power supply is plugged into the BD Veritor™ Plus Analyzer, examine the AC Adapter plugged into the Analyzer. The supplied AC Adapter has a BD label near the coaxial power connector and conforms to the sample labeling in section 2.3.

6.4 Result Messages

The following codes and messages may be displayed on the display window of the BD Veritor™ Plus Analyzer and/or in the RECORD.CSV file that may be downloaded from the BD Veritor™ Plus Analyzer (see section 4.8 for download instructions).

When the system encounters an alert or error condition, the result number and message is usually displayed.
Contact numbers for BD technical services and local BD representatives outside of the United States appear on page 215.

Result Number	Message	Possible Cause(s)	Corrective Action(s)
SYSTEM ALERTS			
01	Varies according to test type	Successful test execution	None required.
02	POSITIVE CONTROL LINE INVALID CONSULT MANUAL	The positive assay control line is outside the normal range.	On first occurrence, repeat the assay. On subsequent occurrences, contact BD*.
03	TEST WAS REMOVED REINSERT TEST DEVICE	The test device was removed while the run was in progress.	That run is invalidated. The device may be re-inserted to begin a new run.
05	BACKGROUND INVALID CONSULT MANUAL ERROR 05	The assay background area is outside of the normal range.	Remove the device and cycle the instrument power. <ul style="list-style-type: none"> If the instrument displays an error during power-on, contact BD*. If no error message appears, obtain an additional sample and repeat the assay. If the issue recurs, contact BD*.
06	ANALYZER FAILURE CALL BD ERROR 06	A hardware anomaly has occurred.	Contact BD*.
07	ANALYZER FAILURE CALL BD ERROR 07	A hardware anomaly has occurred.	Contact BD*.
08	ANALYZER FAILURE CALL BD ERROR 08	A high ambient light condition has been detected or a hardware anomaly has occurred.	On first occurrence, move the unit to a lower ambient light condition. In case of persistent recurrence, contact BD*.
09	ANALYZER FAILURE CALL BD ERROR 09	A hardware anomaly has occurred.	Contact BD*.
11	ANALYZER FAILURE CALL BD ERROR 11	A hardware anomaly has occurred.	Contact BD*.
12	NEGATIVE CONTROL LINE INVALID CONSULT MANUAL	Cross-reacting material present.	On first occurrence, repeat the assay. On subsequent occurrences, contact BD*.
* outside of the United States contact your local BD representative			

Result Number	Message	Possible Cause(s)	Corrective Action(s)
13	ANALYZER FAILURE CALL BD ERROR 13	A hardware anomaly has occurred.	Contact BD*.
14	ANALYZER FAILURE CALL BD ERROR 14	A hardware anomaly has occurred.	Contact BD*.
15	ANALYZER FAILURE CALL BD ERROR 15	A hardware anomaly has occurred.	Contact BD*.
16	ANALYZER FAILURE CALL BD ERROR 16	A hardware anomaly has occurred.	Contact BD*.
17	ANALYZER FAILURE CALL BD ERROR 17	A hardware anomaly has occurred.	Contact BD*.
18	ANALYZER FAILURE CALL BD ERROR 18	A hardware anomaly has occurred.	Contact BD*.
19	BACKGROUND INVALID CONSULT MANUAL ERROR 19	The assay background area is outside of the normal range.	Remove the device and cycle the instrument power. <ul style="list-style-type: none"> • If the instrument displays an error during power-on, contact BD*. • If no error message appears, obtain an additional sample and repeat the assay. If the issue recurs, contact BD*.
20	ANALYZER FAILURE CALL BD ERROR 20	A hardware anomaly has occurred.	Contact BD*.

* outside of the United States contact your local BD representative

Result Number	Message	Possible Cause(s)	Corrective Action(s)
21 / 22	ANALYZER EXPIRED NO DAYS REMAIN	The shelf life of the BD Veritor™ Plus Analyzer has expired. The maximum shelf-life of the BD Veritor™ Plus Analyzer from the date of manufacture is 34 months (regardless of time of first use). or The operational life of the BD Veritor™ Plus Analyzer has expired. The maximum number of months from first use is 24 months.	Contact your distributor to order a new BD Veritor™ Plus Analyzer*.
23	ANALYZER EXPIRED NO TESTS REMAIN	The operational life of the BD Veritor™ Plus Analyzer has expired. The maximum number of tests performed by a BD Veritor™ Plus Analyzer is 10,000 runs.	Contact your distributor to order a new BD Veritor™ Plus Analyzer*. Outside of the United States contact your local BD representative
24	INVALID BARCODE ON TEST DEVICE REPEAT TEST WITH NEW TEST DEVICE	The test device barcode cannot be read correctly.	Repeat the assay. On subsequent occurrences, contact BD*.
25	VERIFY FAIL CONSULT MANUAL ERROR 25	The result of a BD Veritor™ System Verification Cartridge run is unsuccessful.	Ensure that no foreign matter is present on the BD Veritor™ System Verification Cartridge assay area, and repeat the run. On subsequent occurrences, contact BD*.
26	RESULT INVALID FLU A: + FLU B: + CONFIRM RESULT ERROR 26	The BD Veritor™ Plus Analyzer has detected potential interference in the test line area.	Specimens generating a "Result Invalid" with Error 26 should be retested. Upon retesting, if the specimen produces a "Result Invalid" the user may want to consider other methods to determine whether the sample is positive or negative for influenza virus.
27	ANALYZER FAILURE CALL BD ERROR 27	A hardware anomaly has occurred.	Contact BD*.
28	ANALYZER FAILURE CALL BD ERROR 28	A hardware anomaly has occurred.	Contact BD*.
29	ANALYZER FAILURE CALL BD ERROR 29	A hardware anomaly has occurred.	Contact BD*.

* outside of the United States contact your local BD representative

Result Number	Message	Possible Cause(s)	Corrective Action(s)
30	ANALYZER FAILURE CALL BD ERROR 30	A hardware anomaly has occurred.	Contact BD*.
31	UPDATE ERROR 31 TURN OFF POWER REMOVE USB DEVICE	The BD Veritor™ Plus Analyzer update has not been completed successfully.	Repeat the firmware update procedure as described in the instructions included with the update USB flash drive. On subsequent occurrences, contact BD*.
32	UPDATE ERROR 32 TURN OFF POWER REMOVE USB DEVICE	The BD Veritor™ Plus Analyzer update has not been completed successfully.	The firmware update is not compatible with the Analyzer. Contact BD to obtain the correct update version*.
33	UPDATE ERROR 33 TURN OFF POWER REMOVE USB DEVICE	Defective USB flash drive.	Repeat the firmware update procedure as described in the instructions included with the update USB flash drive. On subsequent occurrences, contact BD*.
34	Test name and result followed by: QC FAIL	The BD Veritor™ Plus Analyzer has detected that the QC test procedure has not been completed successfully, because: <ul style="list-style-type: none"> A mismatch was detected between the QC swab and test device (for example, a Strep A QC swab was run with an RSV test device). The test outcome did not match the expected outcome (for example, a Strep A positive result was expected but a Strep A negative result was detected). 	Ensure that the same test device type has been used across all components (for example, a BD Veritor™ System Test Flu A+B test device must be used with the QC swab appropriate for the test device). Ensure that the test device inserted was developed from the QC swab that was scanned with the barcode scanner. On first occurrence, repeat the QC procedure. On subsequent occurrences, contact BD*.
35	CONNECTIVITY MODULE ERROR	A hardware or software failure of the connectivity module.	Contact BD for replacement module*.
none	UPDATE FILE NOT FOUND REMOVE USB DEVICE	Defective USB flash drive or non-supported device inserted into BD Veritor™ Plus Analyzer USB port.	If attempting to update the Analyzer, confirm that the USB flash drive has been provided by BD. Repeat the firmware update procedure as described in the instructions included with the USB flash drive. On subsequent occurrences, contact BD*.
none	INVALID BARCODE REINSERT TEST DEVICE	An invalid kit lot barcode was detected by the BD Veritor™ Plus Analyzer.	Ensure that the barcode scanned during this step is from a valid BD Veritor™ Test Kit box label. For kits with two barcodes, scan only the upper barcode that starts with (17).

* outside of the United States contact your local BD representative


Result Number	Message	Possible Cause(s)	Corrective Action(s)
none	INVALID KIT LOT BARCODE	An invalid kit lot barcode was detected by the BD Veritor™ Plus Analyzer.	Ensure that the barcode scanned during this step is from a valid BD Veritor™ Test Kit box label. For kits with two barcodes, scan only the upper barcode that starts with (17).
none	OPERATOR ID SCAN TIMEOUT REINSERT TEST DEVICE	No barcode was scanned within 30 seconds.	Reinsert test device and scan all barcodes within 30 seconds of the display prompt.
none	SPECIMEN ID SCAN TIMEOUT REINSERT TEST DEVICE	No barcode was scanned within 30 seconds.	Reinsert test device and scan all barcodes within 30 seconds of the display prompt.
none	LOT NUMBER SCAN TIMEOUT REINSERT TEST DEVICE	No barcode was scanned within 30 seconds.	Reinsert test device and scan all barcodes within 30 seconds of the display prompt.
none	XXX DAYS ANALYZER LIFE REMAINING	The shelf or operational life of the BD Veritor™ Plus Analyzer is nearing expiration. The maximum shelf-life of the BD Veritor™ Plus Analyzer from the date of manufacture is 34 months (regardless of time of first use). The maximum number of months from first use is 24 months.	Contact your distributor to order a new BD Veritor™ Plus Analyzer*. The duration of the message is two seconds. To clear the message immediately, insert a test device to begin a new run.
none	XXX READS REMAINING	The BD Veritor™ Plus Analyzer is nearing its usage expiration. The maximum number of tests performed by a BD Veritor™ Plus Analyzer is 10,000 runs.	Contact your distributor to order a new BD Veritor™ Plus Analyzer*. The duration of the message is two seconds. To clear the message immediately, insert a test device to begin a new run.
none	Flashing Plug icon 	Failure of the internal battery. Battery Not Charging: In the event the battery icon does not indicate the battery is charging while the power supply is plugged into the BD Veritor™ Plus Analyzer, examine the AC Adapter plugged into the Analyzer. The supplied AC Adapter has a BD label near the coaxial power connector and conforms to the sample labeling in section 2.3.	Contact BD*. Note: the Analyzer may continue to be used while plugged into a facility power source.

* outside of the United States contact your local BD representative

Result Number	Message	Possible Cause(s)	Corrective Action(s)
none	WALKAWAY MODE REMOVE TEST DEVICE INSERT WHEN PROMPTED	Test device inserted at the incorrect workflow step during Walk Away operation mode.	Current run aborted. Remove test device and immediately restart the Walk Away test sequence. Insert the test device after the display window prompt to do so.
none	TEST DEVICE WAS REMOVED REPEAT ASSAY	Test device was removed before the test run was complete during Walk Away operation mode.	Current run aborted. The test must be repeated with a new test device because the elapsed development time of the assay is unknown.
* outside of the United States contact your local BD representative			

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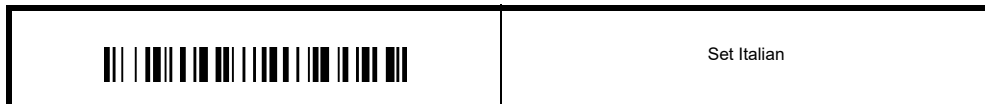Australian Sponsor
Becton Dickinson Pty Ltd.
66 Waterloo Road
Macquarie Park NSW 2113
Australia

New Zealand Sponsor
Becton Dickinson Limited
14B George Bourke Drive
Mt. Wellington Auckland 1060
New Zealand

Symbol Glossary/Glossaire des symboles/Symbolverzeichnis/Legenda dei simboli/ Glosario de símbolos/Symbolförklaring

Manufacturer / Fabricant / Hersteller / Fabricante / Fabricante/ Tillverkare

Use by / Date de péremption / Verwendbar bis / Usare entro / Usar antes de / Använd före

YYYY-MM-DD / YYYY-MM (MM = end of month)
AAAA-MM-JJ / AAAA-MM (MM = fin du mois)
JJJJ-MM-TT / JJJJ-MM (MM = Monatsende)
AAAA-MM-GG / AAAA-MM (MM = fine mese)
AAAA-MM-DD / AAAA-MM (MM = fin del mes)
AAAA-MM-DD / AAAA-MM (MM = slutet av månaden)

REF

Catalog number / Numéro catalogue / Katalognummer / Numero di catalogo / Número de catálogo / Kataloški broj

EC REP

Authorized Representative in the European Community / Représentant autorisé pour la Communauté européenne / Autorisierter Vertreter in der Europäischen Gemeinschaft / Rappresentante autorizzato nella Comunità Europea / Representante autorizado en la Comunidad Europea / Auktoriserad representant i Europeiska gemenskapen

IVD

In Vitro Diagnostic Medical Device / Dispositif médical de diagnostic in vitro / Medizinisches In-vitro-Diagnostikum / Dispositivo medicale per diagnostica in vitro / Dispositivo médico para diagnóstico in vitro / Medicinteknisk produkt för in vitro-diagnostik

This only applies to US: "Caution: Federal Law restricts this device to sale by or on the order of a licensed practitioner." / S'applique uniquement aux États-Unis : "Caution: Federal Law restricts this device to sale by or on the order of a licensed practitioner." / Vale solo per gli Stati Uniti: "Caution: Federal Law restricts this device to sale by or on the order of a licensed practitioner." / Gilt nur für die USA: "Caution: Federal Law restricts this device to sale by or on the order of a licensed practitioner." / Sólo se aplica a los EE.UU.: "Caution: Federal Law restricts this device to sale by or on the order of a licensed practitioner." / Detta gäller enbart USA: "Caution: Federal Law restricts this device to sale by or on the order of a licensed practitioner."

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